

SMART NEBULIZER

• 28123 SMART NEBULIZER

Easy-to-use single button operated compressor nebulizer for home use. Light and compact design.


Reliable pump with service life exceeding 1000 hours.

It ensures highly efficient drug delivery for the treatment of asthma, allergies, and various respiratory diseases.

Supplied with tube, adult and child masks, mouthpiece, 2 air filters and multilingual manual in GB, FR, IT, ES, PT, DE, GR, RO, DK, HU, HR, BG.

• 28059 RESPIRATORY ACCESSORIES KIT - spare

(adult and child mask, 5 filters, mouthpiece, tubing)

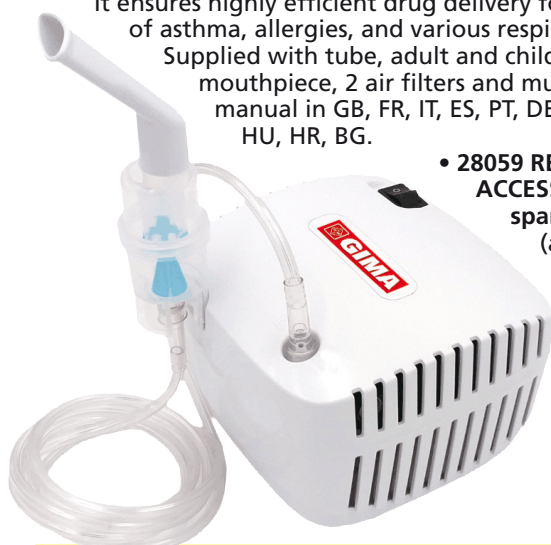

TECHNICAL SPECIFICATIONS

Max. capacity:	10 ml
Air flow:	≥ 7 l/min
Operating mode:	20 min ON, 40 min OFF
Max. pressure:	210-400 Kpa
Nebulizing pressure, rate	60~180 Kpa; approx. 0.35 ml/min
Particle size / FPF:	< 5µm / ≥60%
Altitude operating range:	0-2000 m
Noise level:	approx. 55 dB
Voltage:	AC 230 V, 50 Hz, 0.7A; ≤ 80W
Size - weight:	140x140x100 mm - 1.1 kg
Safety:	Class II, Type BF

MAGIC NEB - MESH NEBULIZER

• 28065 MAGIC NEB WEARABLE MESH NEBULIZER

Pocket size and wearable ultrasonic nebulizer for treatment of respiratory diseases. Reliable quality, soundless operation and fine particle size. Self-cleaning mode: at the end of the treatment, the device mists pure water cleaning the mesh from medication residues.

- battery level indicator

- type-C charge port

- anti-corrosion material

- one-button operation

Supplied with adult mask,

child mask, inhalation

mouthpiece, built-in

battery, type C cable and

multilingual manual

GB, FR, IT, ES, PT,

DE, DK, RO,

HR, GR.



STANDARD ACCESSORIES



CE 0051

28065

TECHNICAL SPECIFICATIONS

Ultrasonic frequency:	approx. 110 kHz	Head life:	180 hours
Nebulizing rate:	≥0.2 ml/min	FPF:	≥ 60%
MMAD:	approx 3.5 µm		
Max./min. liquid volume:	6 ml / 2 ml		
Shutdown:	automatic without liquid		
Power supply:	built-in lithium battery 3.7 V, 300 mAh		
Operating time:	> 60 min.		
Battery life:	500 cycles		
Size - weight:	54x36x44 mm - 45 g		
Liquid protection:	IPX5		

CLINEA NEBULIZER

• 28058 CLINEA NEBULIZER - piston

High-end, low-noise compressor nebulizer designed to provide efficient and quiet aerosol therapy.

Ideal for professional healthcare settings.

Featuring innovative brushless DC motor technology, this device combines power, reliability and long-lasting performance to deliver a comfortable and effective treatment experience.

Equipped with an adjustable pressure function, the nebulizer allows users to adjust the airflow intensity according to their specific needs or medical recommendations.

Supplied with adult mask, child mask, bulb, inhalation mouthpiece, inhalation tube, air filter, adapter and multilingual manual GB, IT




28058

• 28059 RESPIRATORY ACCESSORY KIT (adult and child mask, 5 filters, mouthpiece, tubing) - spare

TECHNICAL SPECIFICATIONS

Operating mode:	continuous use	FPF:	≥60%
Nebulizing flow:	≥4 L/min	Nebulizing Rate:	≥0.2 ml/min
MMAD:	<5 µm	Max capacity:	10 ml
Noise level:	≤55 dB	Operating time:	> 4 hours
Expected pump life:	1000 hours	Size:	168.5x116x99 mm
Voltage:	100~240V 50~60Hz, 700mA Max		
Input:	DC 12V 1.5 A		

REMEDY RHINITIS RELIEVER

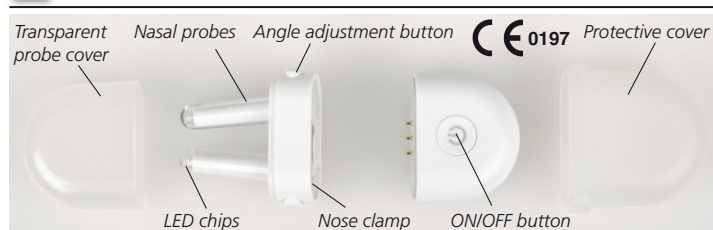

• 28107 REMEDY RHINITIS RELIEVER

This device offers an innovative and effective treatment for the relief of symptoms associated with allergic rhinitis like watery eyes, runny nose and sneezing after a few treatments.

Safe, quick, portable and easy to use. It works by using phototherapy of two different wavelengths: one (infrared light) produces small amounts of heat to improve blood flow and circulation and the second wavelength (red light) suppresses the cells that release histamine and lowers inflammation, thereby alleviating the irritating symptoms of allergic rhinitis.

- 3 minute treatment

- wireless nasal probes

- effective relief of allergy symptoms

- red light therapy

- one button operation

Medical Device Class IIa

Supplied with nose clamp,

2xCR 2032 batteries, protective cover

and multilingual manual GB, FR, IT, ES,

PT, DE, NL, SE, RO, GR.



TECHNICAL SPECIFICATIONS

Treatment system:	dual wavelength phototherapy
Wavelength:	652 nm (red light); 940 nm (infrared light)
Operating time:	3 minutes
Light power:	3-6 mW per nostril
Power supply:	2xCR2032 alkaline batteries, 3V
Size - weight:	70x39x17 mm - 21 g (without batteries)
Classification:	internally powered, Type BF applied part, IP22

