

HEINE DELTA 30 UV DERMATOSCOPE

• 31144 HEINE DELTA 30 UV DERMATOSCOPE - K-240.28.305

HEINE Delta 30 UV Dermatoscope unites ultraviolet, polarised, and non-polarised illumination in one advanced dermatoscope.

Twelve UV LEDs with a wavelength of 365 nm provide uniform illumination of the examined area.

Specialized filters ensure pure, consistent UV light while eliminating stray wavelengths, enhancing fluorescence contrast for clear and reliable visualization.

- the angled LEDHQ arrangement enhances visual depth, optimises clarity, and reduces reflections

- 32 mm lens delivering a "real" 30 mm field of view

- 10x magnification for detailed examination

- natural colour rendering for enhanced contrast and accurate assessment

- integrated dioptr adjustment for more precise focusing

Its intuitive handling and seamless integration of multiple illumination options streamlines examinations and supports confident clinical decision-making.

Supplied with contact plate with scale, rechargeable Li-ion battery, USB-C power supply with cable, 4 adaptors and case.

• 31147 HEINE CHARGING STATION - X-002.99.212 - optional

• 31148 HEINE SMALL CONTACT PLATE 8 mm - K-000.34.206 - optional

• 31149 UNIVERSAL SMARTPHONE ADAPTOR MAG - K-000.34.273 - optional

• 31151 HEINE STANDARD CONTACT PLATE 30 mm - K-000.34.220 - spare

Three illumination options in one device
Toggle easily between UV, polarised, and non-polarised illumination

5 YEARS WARRANTY

Smooth and easy-to-clean surface

31144
Front

Achromatic HEINE optical system for an optimal visual field

Long-lasting and robust for high reliability
- Ultra-light aluminium frame for strength and durability
- LEDHQ with virtually unlimited lifetime and consistently strong performance

Custom light intensity
- Side-mounted button with three degrees of luminous intensity
- Adjustable for both white light and UV illumination

MADE IN GERMANY

UV illumination for fluorescence visualization

Smart power source options

Intuitive and ergonomic

31144
Back

TECHNICAL SPECIFICATIONS	
Dioptr:	-4 to +4 dpt
Illuminance:	22.000 lux in 15 mm distance without contact plate
Colour temperature:	5700 K
Color rendering index (CRI):	92
Working distance:	15 mm in non-contact mode or contact to skin
Power supply:	Li-ion built-in rechargeable battery
Input:	USB-C, 5 V DC, 1.2 A or charging station
Charging time:	150 min
Size - weight	195x55x70 mm - 280 g
Protection:	IP 20

HEINE K180 OTO-OPHTHALMO

• 31739 HEINE K180 - LED F.O. OTO-OPHTHALMO SET - 2.5 V - A-279.10.118

Otoscope and ophthalmoscope K180 with BETA battery 2.5 V handle and a set of single use specula. Supplied in hard case. Housing made of lightweight, shockproof polycarbonate.

• 31517 HEINE SANALON S REUSABLE EAR SPECULUM for K180, K200 - set of 4 (2.4/3/4/5) - B.000.11.111

• 31742 HEINE K180 LED OPHTHALMOSCOPE - 2.5 V

- C-182.10.118

27 lenses from -35D to +40D

5 different apertures in 2 variants.

31742



5 YEARS WARRANTY



MADE IN GERMANY

31739

HEINE BETA 200 F.O. OTOSCOPES



MADE IN GERMANY


31724



34449



See page 376

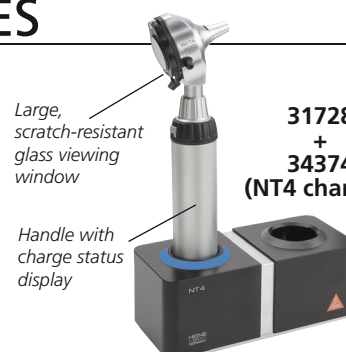

Large, scratch-resistant glass viewing window

Handle with charge status display

31728 + 34374 (NT4 charger)

Distal fiber optic (F.O.) illumination with LED technology. Airtight window seal for reliable test of tympanic mobility

• 31724 HEINE BETA 200 LEDHQ F.O. OTOSCOPE - 2.5V - B-141.10.118

Exceptional durability combined with optimum 3X magnification.

Particularly bright illumination thanks to LEDHQ.

Solid metal housing.

Supplied in hard case with BETA battery 2.5V handle with rheostat, a set of single use specula and a set of reusable specula.

• 34449 HEINE NT 300 CHARGER

• 34448 HEINE RECHARGEABLE BATTERY - X-007.99.383



31729

• 31729 HEINE BETA 200 LEDHQ OPHTHALMOSCOPE - 2.5 V - C-144.10.118



6 standard diaphragms: slit, fixation star with polar co-ordinates, cobalt blue filter, large spot, small spot, hemi spot



6 diaphragms with green filter
+ in 1 D steps 1 - 10 | 15 | 20 | 40 | D
- in 1 D steps 1 - 10 | 15 | 20 | 25 | 35 | D