

SPINAL IMMOBILIZERS

- **34032 SPINAL IMMOBILIZER - adult M (150-175 cm height)**
- **34062 SPINAL IMMOBILIZER - adult L (170-185 cm height)**

Designed to secure the patient, the surrounding is suitable for individuals with rib injuries.

Strips, available in different colours, can be quickly fixed. It is equipped with nylon buckles and safety belts.

It is recommended that when used in conjunction with a neck fixation, the patient should be in sitting position, before being transferred to a PE stretcher.


Supplied with forehead strap for head fixation and safety strap for chin fixation. Max. load capacity: 175 kg.

Weight: 3 kg. Unfolded size: 79x82 cm.

Folded size: 84x25x9 cm (34032)

Weight: 4.5 kg. Unfolded size: 88x86 cm.

Folded size: 88x28x9 cm (34062)

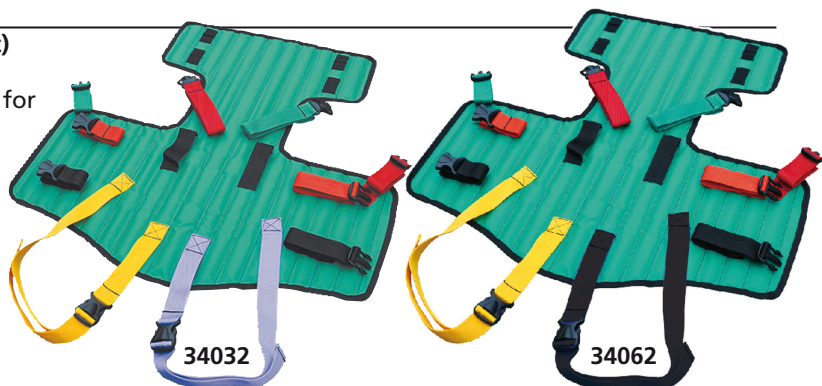

Multilingual manual: GB, FR, IT, ES, PT, DE, NL, DK, SE, FI, LT, LV, EE, PL, CZ, SK, HU, SI, HR, RO, BG, GR, Arabic


TRACTION SPLINT

- **55990 TRACTION SPLINT - adult**

Traction splint relieves tissue pressure on fractured leg by avoiding further haemorrhages, shock and secondary damages to nerves, bones, muscles and blood vessels. Adjustable length and automatic traction mechanism.

Made of aluminium alloy and nylon leg straps. Supplied with carrying bag.

Size: 134x19x21 cm. Folded size: 93x20x7 cm. Weight: 2 kg.

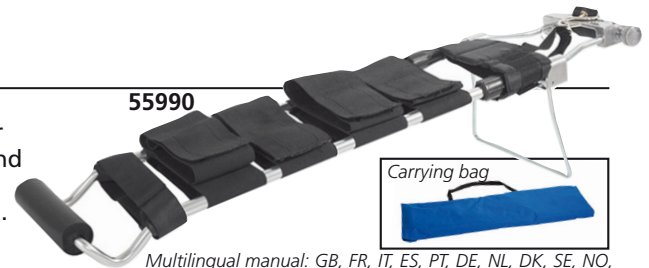

Multilingual manual: GB, FR, IT, ES, PT, DE, NL, DK, SE, NO, FI, LT, LV, EE, PL, CZ, SK, HU, SI, HR, RO, BG, GR, Arabic

VACUUM SPLINT KIT

- **34690 VACUUM SPLINT KIT**

The kit includes vacuum splints for arm, elbow and leg, nylon bag and pump. Made of high quality and durable TPU material, filled with harmless EPS particulate material. Suitable for X-ray and CT scan. The firmness of the splints can be adjusted by pumping air according to the situation.

Quick and easy softness and hardness adjustment of the splint. Packaging size: 50x28x22 cm



Carrying bag and 0.25 l pump included



GIMA CERVICAL COLLARS



- **34615 GIMA CERVICAL COLLAR - pediatric**
- **34616 GIMA CERVICAL COLLAR - adult**

These cervical collars are designed for maximum versatility and patient safety. Adjustable size make it easy to customize the fit for most patients. They feature a functional tracheal aperture: this port provides immediate access for treatment of throat and trachea. Fully compatible with X-ray, CT, and MRI scans. Non-toxic, latex-free and hypoallergenic components.

Size: 45x14x1 cm (34615) 56x19x1 cm (34616)

Weight: 0.1 kg. Multilingual instructions: GB, FR, IT

ROLLED SPLINTS



- **55900 ROLLED SPLINT 92x11 cm**
- **55901 ROLLED SPLINT 100x11 cm**

This lightweight, waterproof splint is easily mouldable into structural arcs to stabilize fractures, thanks to its foam and aluminium core. It can be cut to any length and only requires tape to secure it. The closed-cell padding ensures patient comfort. Its durable design allows it to be washed and reused for pre- and post-cast care. X-ray translucent. Thickness: 0.5 mm. Multilingual instructions