

Supplied with box and manual in:
GB, FR, IT, DE, ES, NL, PL.
From web you can download:
GR, PT, TH, Arab.

HISTOFREEZER® CRYOSURGICAL SYSTEM 0543

- 30587 HISTOFREEZER® - kit with applicator 2 mm 170 ml aerosol bottle and 60 applicators 2 mm.
- 30588 HISTOFREEZER® - kit with applicator 5 mm 170 ml aerosol bottle and 50 applicators 5 mm.
- 30586 HISTOFREEZER® - large kit mixed applicators 170 ml aerosol bottle, 30 applicators 2 mm, 30 applicators 5 mm.
- 30585 HISTOFREEZER® - mini kit mixed applicators One 80 ml aerosol bottle, 15 applicators 2 mm, 15 applicators 5 mm.


The Histofreezer® portable cryosurgical system is a low cost, highly effective cryosurgical agent. This therapeutic treatment makes use of local freezing for the controlled destruction or removal of living tissue.

Effective use on 9 benign lesion types


- Verruca Vulgaris (Common Warts)
- Verruca Plantaris (Plantar Warts)
- Molluscum Contagiosum
- Lentigo (Age Spots)
- Condyloma(Genital Warts)
- Seborrhic Keratosis
- Acrochordon (Skin Tags)
- Verruca Plana (Flat Warts)
- Actinic keratosis (Facial and Non-Facial)




2 mm disposable applicator for small lesions
5 mm disposable applicator for large lesions

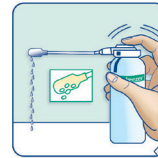

CHARGE
custom applicator
press valve until
droplets start to form.

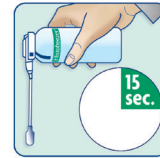

INVERT
canister
for 15 seconds for optimal
freezing temperature.

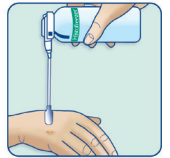

APPLY
a frozen
applicator to lesion
for 15-40 seconds.

Saturated 2 mm applicators will stay cold for 90 seconds.
Saturated 5 mm applicators will stay cold for 120 seconds.

CHANGE-A-TIP™ CAUTERIES WITH HIGH TEMPERATURES


BOVIE CAUTERIES - made in USA

- interchangeable sterile tips
- replaceable batteries
- reusable handles
- no sterilization necessary



30590
(Fine tip, not sterile, included)

| GIMA code | INTERCHANGEABLE TIP CAUTERIES MINIMUM ORDER: 10 TIPS | Recommended for: |
|-----------|---|------------------|
| 30590 | Cautery handle "AA" (tip included) | All tips |
| 30593 | High temp loop tip - sterile | DER/ENT |
| 30594 | High temp 2" flex loop tip - sterile | GYN/ENT |
| 30595 | High temp 5" loop tip - long size - sterile | GYN/AECOLOGY |
| 30596 | High temp fine tip - sterile | DER/CARDIO |

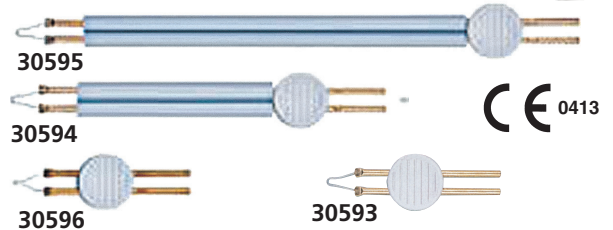


 0413

REUSABLE HIGH TEMPERATURE CAUTERY WITH REPLACEABLE TIP

- CAUTERIES - 1,200°C (2,200°F) - made in Italy
- interchangeable tips with nickel-chrome wire
- replaceable batteries
- reusable handles

Wide application in gynaecology, orthopedy, ear/nose/throat, dermatology, plastic and reconstruction surgery.

| GIMA code | INTERCHANGEABLE TIP CAUTERIES MINIMUM ORDER: 10 TIPS |
|-----------|---|
| 30598 | Cautery handle "AA" (tip not included) |
| 30599 | High temp long fine tip - sterile - box of 10 |
| 30600 | High temp large tip - sterile - box of 10 |
| 30601 | High temp fine tip - sterile - box of 10 |


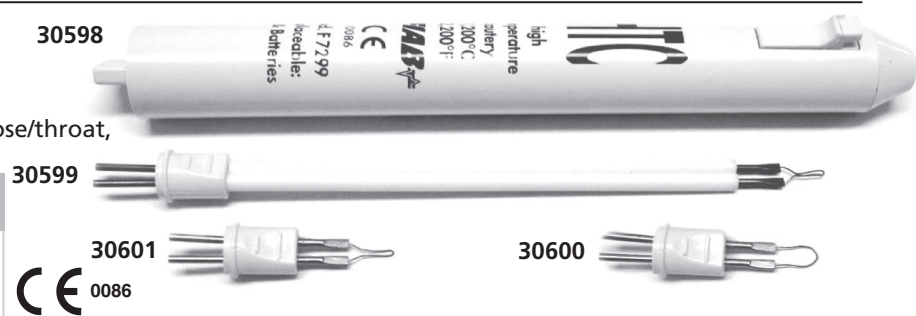


 0086

HIGH TEMPERATURE ELECTROCAUTERIES

- 30580 ELECTROCAUTERY - 1,200 °C - fine tip - made in Italy
Electrocautery is a sterile bipolar battery operated disposable device for the cauterization of small vessels; the temperature of the tip is 1,200°C. With a fine tip required in most of the surgeries.
- 30582 GYNAECOLOGY CAUTERY - 1,200°C - fine tip
Same as 30580 but with the tip fixed on a 12 cm long shaft.
- 30581 ELECTROCAUTERY - 600 °C - fine tip
Same as 30580 but with lower temperature = 600°C. Especially designed for ophthalmology.
- 30583 EMERGENCY ELECTROCAUTERY - 800°C - Thick tip
Same as 30580 but with "thick tip" for the perforation of nails in case of haematoma.

 0086

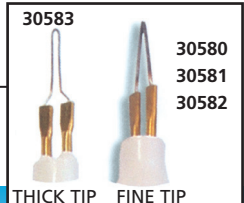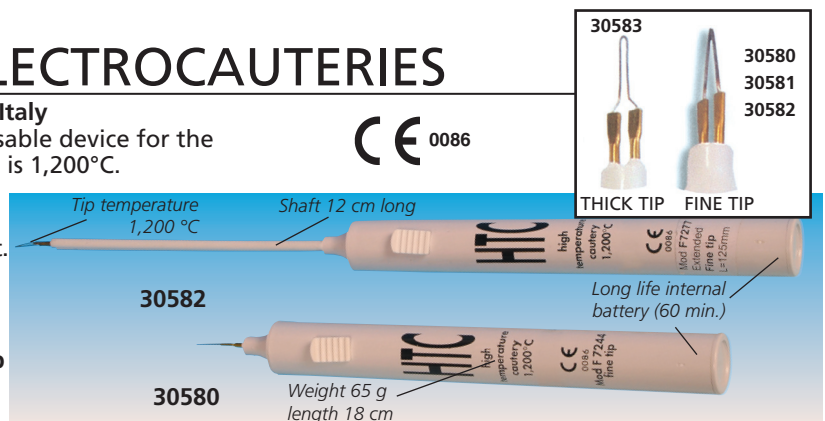

THICK TIP FINE TIP

30582

30580

Weight 65 g
length 18 cm

Long life internal battery (60 min.)